

Complete  
and versatile

Reprocessing can be as  
simple and efficient: with  
the optimum equipment  
accessories.

1 – Upper  
basket

Flexible for  
inserts  
and height  
adjustable.



2 – Insert

Intended for  
different  
instruments  
and for  
flexible use.



3 – Lower  
basket

For holding  
mesh trays  
as well  
as various  
inserts.



4 – Insert

For  
accommodating  
various  
instruments and  
insertable for  
upper and lower  
basket.




»starter pack« 1–6

includes the perfect range  
of basic equipment.

»injection pack« 3–7

complemented by an  
additional equipment perfectly and  
particularly suitable for handpieces.

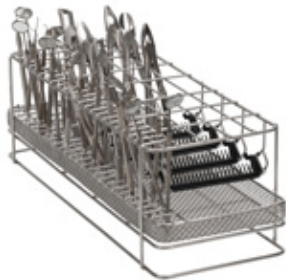

5 – Insert

Different  
compartmentds  
for optimum  
loading  
of upright  
instruments.



6 – Insert

Equipped  
with different  
compartments  
for optimum  
loading.



7 – Upper  
basket with  
injector

For various  
instruments with  
injector nozzles.

Technical  
data

CE  
0297

Type:	Teon	Teon +
Power supply:	400 V AC; 50–60 Hz/230 V AC; 50–60 Hz	
Current:	16 A	
Overall dimension (w × h × d):	600 × 835 × 600 mm	
Space required for installation (w × h × d):	598 × 835 × 1,200 mm	
Weight (empty):	74 kg	78 kg
Smart Dry:	•	–
Dry Plus:	–	•
HEPA filter:	–	•
Cold water connection:	•	•
Hot water connection:	–	•
Demineralized water connection:	•	•
Auto-close – automatic door lock:	•	•
No. of integrated dispenser pumps for liquid products:	2	1
In-door dosing device for rinse:	•	•
Salt tank – front door access:	•	•
Number of program:	5	6
No. of customizable program:	2	2
Cycle time*:	49'	74' including active air-drying
* time for Vario TD Dental program		

The thermo washer disinfectors Teon and Teon + were designed, certified and validated with the most stringent directives and standards:

93/42/EEC Medical Device Directive	EN 61010-2-040 Safety requirements	DIN EN ISO 15883-1/2 Washer-Disinfector
EN 61326-1 Electromagnetic compatibility	–	EN1717 Protection against pollution of potable water

W&H Sterilization Srl.  
Via Bolgara, 2  
24060 Brusaporto (BG), Italy  
t +39 035-66 63 000  
f +39 035-50 96 988  
office.sterilization@wh.com  
wh.com

#whdentalwerk  
f @ in ▶  
video.wh.com

Now at  
your dealer  
or **wh.com**



Top  
performer  
for high cleaning  
standards



New!

teon  
teon +

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20492 AEN Rev. 001 / 04.03.2019  
Subject to alterations



## Teon

### Next step for reprocessing

The new Teon is a highly efficient thermal washer disinfecter. Thanks to its dedicated programs it offers the best cleaning and disinfection for your dental practice. Teon improves the job organization of your hygiene room, reduces the reprocessing time and increases the protection against cross contamination.



Make it the way you need it.  
For example with the  
»starter pack« inside.

teon

## High standard for cleaning and reprocessing



## Teon +

### The best approach for reprocessing

Teon +, the premium W&H thermal washer disinfecter, ensures the most demanding clinics an excellent automatic cleaning system for instruments, including handpieces. Teon + offers full traceability to meet current regulatory requirements. Its patented innovative active air-drying system offers a new standard in terms of quality and performance.



Complete your device with the  
equipment you really need,  
such as with the »injection pack«.

teon +